gorenje

Elektronska osebna tehtnica
Elektronska osobna vaga
Elektronska osobna vaga
Electronic personal scale
Cântar de bucătărie electronic
Электронни лични кантар
Електронні ваги
Elektroniczna waga osobista
Elektronická osobná váha
Személymérleg
Электронные весы
Bilancia pesapersona

OT 180 GB
OT 180 GW

Navodilo za uporabo SI

Upute za uporabu BIH HR

Uputstva za upotrebu BIH SRB MNE

Instructions for Use GB

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gorenje

ATTENTION

Thank you for choosing Healthcare OT 180 GB-GW scale. To ensure a proper operation procedure, please carefully read through the instruction manual.

This appliance is marked according to the European directive 2002/96/EC on Waste Electrical and Electronic Equipment (WEEE). This guideline is the frame of a European-wide validity of return and recycling on Waste Electrical and Electronic Equipment.

HELPFUL TIPS

Your scale is a precision electronic instrument. One of the hardest things to weigh is the human body because it is always in motion. To ensure reliable readings, always try to place your feet on the scale in the same place and stand still.

Treat your scale with care. Do not drop it or jump on it. The scale is durable, but misuse or abuse may render the electronic sensors inoperable.

Never submerge the scale under water. To clean the scale surface, use either a damp cloth or glass cleaner. Do not open or take apart the scale or its components.

HOW TO OPERATE

1. Open the battery cover at the backside of the scale.
2. If there are batteries in the case, take the insulator out. Otherwise insert 1x CR 2032 batteries to the battery case. Then close the battery cover.

Figure 1

Figure 2


Figure 3

3. Select the weight unit you want among kg/lb/st.

Kg/lb/St

4. Place the scale on a hard and even floor.
5. If there is a number other than 0.0 kg displayed, please wait for about 30 seconds until it automatically turns off; if there is nothing on the screen, please overlook next step .
6. With your foot, firmly tap the centre of the scale until "8" displayed .

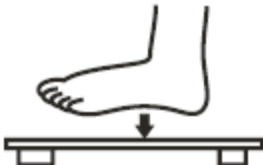7. **Important!** Wait for "0.0 kg" to be displayed before stepping on the scale .

8. Step on the scale (the display will flicker) and your weight measurement will be shown on the screen. The scale will automatically shut off after a few seconds.

TROUBLESHOOTING GUIDE

If you experience any difficulty in use of the scale:

1. Check the battery is correctly fitted.
2. Check that you have selected your preference of kg/lb/st.
3. Check whether the scale is on a flat, level floor and not touching against a wall.
4. Ensure whether you do not step on the scale until the display shows "0.0 kg".
5. If, when you tap the scale to activate it, there is no display at all or "Lo" is displayed- try two new batteries.

Replace battery

6. If "Err" is displayed, the scale has been overloaded.

Weight overload

SPECIFICATION

Weight range Min 5 kg
 Max 180 kg

Graduation 0.1 kg

CAUTION

Keep the scale off water

Keep the scale with care

Don't put it on uneven floor while weighing

Don't step on the edge of the scale

Environment

Do not throw away the appliance with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment.

gorenje

WISHES YOU A LOT OF PLEASURE WHILE USING YOUR APPLIANCE

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