ARDESTO

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IMPORTANT SAFEGUARDS

Please read the safety instructions before using this product. The safety instructions provided here help you use this product safely and correctly, and also to prevent you from injuring yourself or others. When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

READ ALL INSTRUCTIONS BEFORE USING

DANGER-WATER

To reduce the risk of burns, electrocution, fire, or injury to persons:

- Do not place or store the product where it can fall or be pulled into bathtub, sink or toilet.
- Do not place or drop the charger into water or other liquid.
- Do not reach for the charger base that has fallen into water or other liquid. Unplug immediately.
- Do not use while bathing.

DANGER-VOLTAGE

- Do not plug the charger base into an outlet with a voltage other than specified on the bottom of the charger (100-240V).
- Voltage converters do not guarantee compatibility.
- Do not force the plug into an outlet.
- If the plug does not fit easily into the outlet, discontinue use.

MEDICAL-WARNINGS

- If you have had oral surgery or gum surgery within the past two months, consult your dentist before using the product.
 - Consult your dentist if excessive bleeding occurs while using this product or minor bleeding continues to occur after one week use.
- Consult your physician or implant device manufacturer prior to using this product if you have questions concerning the combined use of this product with a pacemaker or other implanted medical device(s).

WARNINGS

To reduce the risk of burns, electrocution, fire, or injury to persons:

- This device can be used by children over 8 years of age and persons with physical or sensory or mental disabilities or with insufficient experience and knowledge, if they are under constant supervision or instructed in the safe use of the device and they understand the possible dangers. Cleaning

and maintenance by the consumer should not be carried out by children without supervision.

- Children should be supervised to ensure they do not play with the product.
- Use this product only for its intended use as described in this manual. This product is designed to clean your teeth, gums and tongue only. Do not use any attachments or brush heads other than those recommended by the manufacturer.
- If your toothpaste includes peroxide, baking soda or other bicarbonate (common in whitening toothpastes), please thoroughly clean your brush head and handle to avoid cracking the plastic.
- Certain liquid mouthwash may damage the device by causing the housing or tank to crack, when using mouthwash, rinse the tank and main unit with water immediately after use. Be especially careful when washing the hose connections at the main unit and the tank.
- Never fill the tank with hot water! The gum lines can be very sensitive and can be easily damaged by hot water.
- Use the unit within an ambient temperature range from zero to 40 ° C.
- Never operate this product if it has a damaged cord, plug or charger. Do not use this product if it is not working properly, if it has been dropped or damaged or dropped into water. Immediately discontinue use if the product appears damaged in any way.
- Keep the product away from heated surfaces.
- If you experience discomfort or pain in your mouth, discontinue use of this product and contact your dentist.
- Do not clean the product in the dishwasher.
- Consult your physician prior to using the toothbrush if you have any medical concerns.
- The toothbrush has been tested and is certified as per ETL regulations.
- Never drop or insert any object into any opening.
- For household use only. Do not use outdoors.
- Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.

SAVE THESE INSTRUCTIONS

EQUIPMENT:

Irrigator - 1

Classic irrigator tip - 2


Irrigator tongue cleaner - 1

Electric Toothbrush - 1

Toothbrush head - 2

Adapter - 1

Instruction manual (contains a warranty card) - 1

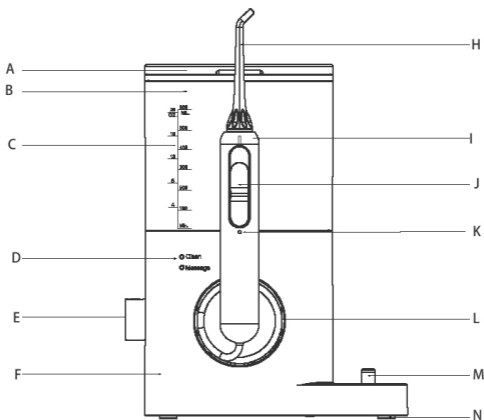ORAL IRRIGATOR**A:** Water Tank Lid**B:** Water Tank**C:** Volume Mark**D:** Mode indicator**E:** Button: Power/Mode**F:** Housing Top**G:** Tank Valve**H:** Nozzle**I:** Release Button**J:** Pause Switch**K:** On/off Mark**L:** Pipe**M:** Charger Base**N:** Rubber Mat**O:** DC Jack Plug**P:** Housing Bottom**G:****O:****P:**

HOW TO USE

GETTING STARTED

Place the device on a flat surface without any sharp objects and make sure the cord is dry and has no damage. The device is powered by a special adaptor with DC jack connector. (Input:100-240V~,50/60Hz;Output:12V,1.5A).

1. Insert the plug into the DC jack on the bottom of base and plug the adaptor cord into an electrical outlet.(Fig.1)
2. Take the Nozzle Case off.
3. Take the Water Tank off to fill with proper liquid, no warmer than 40 °C, do not pour in salt, tooth-paste, or chemicals, then replace the tank onto the top of main unit and press down firmly. (Fig.2-1&Fig2-2)

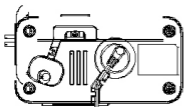
Fig.1

Fig.2-1

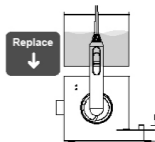
Fig.2-2

INSERTING & REMOVING NOZZLE

1. Inserting Nozzle

Insert tip into the center of the handle, audibly indicating when the nozzle attached to the handle in place. And also try to pull the tip out gently to confirm the tip is locked firmly.(Fig.3)

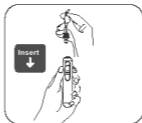
Fig.3

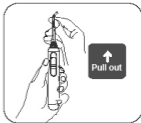
Fig.4

2. Removing Nozzle

Press the Release Button on the top end of the handle, then pull the nozzle out lightly by hand.(Fig.4)

HOW TO USE THE WATER FLOSSER

Do not power on the device before placing the nozzle in your mouth. First use, we suggest you choose the lower pressure, gradually increase the water pressure over time to the setting you need.

1. Hold the handle and point the nozzle to sink, then press down the Main Button to power the device on. (Fig.6)

Press Button again to select the mode you desired. (Clean or Massage), default Clean mode for the first time.

Clean: for cleaning 30 Seconds interval time.

Massage: for cleaning and massaging the gums.

Mode indicator:

White light is on when switching mode;

The Clean mode lamp is above and the Massage lamp is below. (Fig.7-1)

2. Rotate the knob to adjust water pressure, after

setting, push the Pause Switch to "o" position to pause the water flow. (Fig.7-1&Fig.7-2)3. Adjust the direction of nozzle by rotating the nozzle knob;

4. Bend over to the sink and place the tip in your mouth, and then push the pause Switch to «I» position to enjoy irrigating. (Slightly close your mouth to avoid the water splashing.) (Fig.8)

5. Position the water jet flow at 90-degree angle to the teeth and gums. (Fig.9)

6. Move the water jet slowly along the teeth. The unit will continue to work till 2-minute cycle is complete or you turn it off .

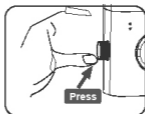
Fig.6

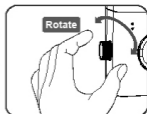
Fig.7-1

Fig.7-2

Fig.8

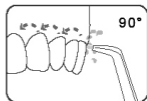
Fig.9

HOW TO END USE

You can end irrigating through the below methods:

1. Auto shut-off turns off after 2- minute cycle is complete;
2. Turn it off by pressing the Main Button;
3. Unplug the adaptor directly when in any emergency.

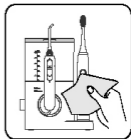
Fig.10

AFTER USE

1. Take off the water tank and dispose of water in it.
2. Push the Pause Switch to " I " position, and power on to drain the water inside the device, then push the Pause Switch to " 0 " position again .
3. Turn off the device and wipe it with a Storage.
4. Soft dry cloth.(Fig.10)
5. Press down the Release Button, and remove the nozzle.
6. Unplug the adapter.

FUNCTIONAL TIP

The number and type of tips depend on delivery set (equipment).

Tongue Cleaner

When you want to clean your tongue, please gently place the tongue cleaner onto the center midline of your tongue surface about halfway back with light water pressure setting first, and gradually increase the pressure as needed. (Fig.11)

Fig.11

Orthodontic Tip

The orthodontic tip is specially designed for people wearing braces: it will clean the hard to reach areas around and in braces.

- 1: Aim water stream at 90 degree to the tooth at the gum-line , start on the last molar (back tooth);
- 2: Slowly slide the tip along with the gum-line , and pausing briefly to use the brush bristles to clean the area between teeth and all around orthodontic braces.
- 3: Repeat until you have finished cleaning (Fig. 12).

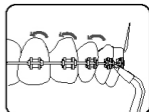
Fig.12

Fig.13

Pocket Tip

The pocket tip is specifically designed for the cleaning of periodontal pockets.

To use the pocket tip, it is a **MUST** to set the device to the lowest pressure setting.

1. Place the soft tip against a tooth at 45- degree angle and gently place the tip under the gum-line, then insert into the pocket. (Fig.13)
2. Turn the unit on and continue tracing along the gum-line until all periodontal pockets are treated. (Fig.14)
3. This tip must be used only with the lowest water pressure to extend life of soft tip as it is easily broken.

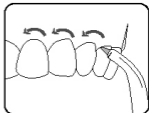
Fig.14

This model only contain classic nozzle and tongue cleaner, other accessories need to be purchased separately.

MAINTENANCE

Use only water or neutral cleaning agent to clean the device. Never use caustic or abrasive agents (e.g., vinegar or lime scale removers), since they may damage the device. Do not wash using hot water of over 40 °C. Keep the device away from high temperature area and avoid direct sunlight.

Main unit

Remove the water tank and drain the remains water inside the device. Clean by using a soft cloth with liquid soap and dry with a soft dry cloth. Do not submerge the main unit in water to clean.

Water tank

Wash with water and wipe with a soft dry cloth. Make sure the water tank is inverted and firmly fixed to store.

Handle

Wash with water and dry up with a soft cloth. Do not bend, pull or twist the water hose.

Nozzle

Wash with water and wipe with a soft dry cloth. Replace once every 6 months for classic nozzle and tongue cleaner, and once every 3 months for other tips.

Charging adaptor

Wipe the metal part of the power plug with a dry cloth once every 6 months.

TROUBLESHOOTING


Problem	Possible Cause	Solution
Unit does not work	Motor not running	Make sure the adapter is firmly connected to the outlet and the unit.
		Make sure the outlet and unit power is on
Water pressure is very low	Dirt stuck in the nozzle	Clean or replace with a new nozzle
	The nozzle is deformed	Replace with a new nozzle
	Dirt stuck in the tank valve	Clean the valve and fill the tank with pure water
Water does not come out	Dirt stuck in the tank valve	Fill the tank up
	Water pressure under the lowest setting	Slide the switch to adjust pressure
	Tank valve not open	Press down the tank, or reassemble the tank valve
Tank leaks	The valve was lost	Replace with a new tank valve
	The valve isn't install in place	Take off the valve and reinstall again
Tear in Pocket tip head	Tip head is worn	Replace with a new nozzle

Simplified Declaration of Conformity:

Shenzhen Risun Technology Co., Ltd. hereby declares that the type of radio equipment (irrigator and electric toothbrush) complies with the Technical Regulations for Radio Equipment. The full text of the declaration of conformity is available on the website at: <https://erc.ua/uk/declarations/>

frequency band(s) in which the radio equipment operates: 148.5 - 5000 kHz
 maximum radio-frequency power transmitted in the frequency band(s) in which the radio equipment operates: minus 15 dB μ A/m
 Software version: software does not apply



TOOTHBRUSH

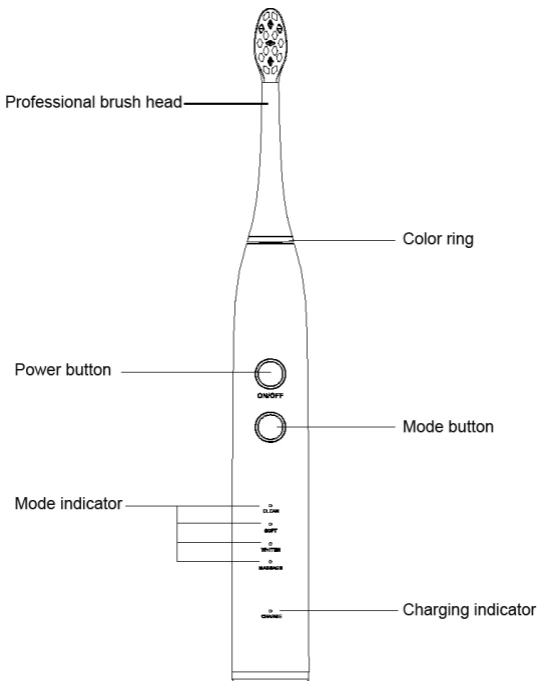GETTING STARTED

Charging the toothbrush

Be sure to turn off the toothbrush before charging. plug the power plug into the outlet. Place the handle in the charger. The charger indicator light will be lit, indicating the handle is charging. After the battery is fully charged, charge indicator light will turn off, indicating the toothbrush is.

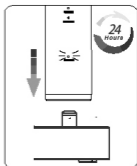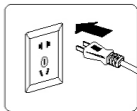Notes:

1. It takes at least 24 hours to charge the battery fully.
2. The toothbrush must be charged for 24 hours, if:
 - 1) this is the first time using the toothbrush;
 - 2) the toothbrush has not been used for more than three months;
 - 3) To keep the battery fully charged, we recommend that you keep the toothbrush in the charger when not in use.

0 During charging, the toothbrush will get warm, but this is normal.

4) The operating duration will become shorter as the battery gets older. If the toothbrush operates only for 1 short time after being fully charged, this means that the battery is nearing the end of its life cycle. Contact service center for support.

HOW TO USE

Slide the brush head onto the handle. Push brush head down until it locks into place on handle.

Note:

1. There is small gap between color ring and the handle. (Fig 1)
2. Wet the bristles and apply a small amount of toothpaste.
3. Place the bristles against the teeth at a slight angle towards the gum line. (Fig 2)
4. Press the mode button to select mode. If you know the mode used last time, you need not select mode again.

There are 4 modes.

CLEAN: Standard mode for superior teeth cleaning.

SOFT: Gently caress, carefully clean sensitive mouth.

WHITEN: Teeth cleaning with polishing and whitening.

MASSAGE: Pulsates for gentle gum stimulation.

Notes:

During brushing, if press Mode button, toothbrush will work at next mode and start another 2 minutes cycle and the mode indicator will light up.

5. Apply light pressure and gently move the brush head slowly across the teeth in a small back and forth motion.

6. In order to maximize cleaning effectiveness, divide your mouth in 4 sections. Start with Section 1 and brush each section for 30 seconds. The brushing action will pause between each 30 second interval, indicating when to move on to the next section.

7. The toothbrush will automatically stop after 2 minutes of use. (Fig3)

8. Pull brush head out from the toothbrush.(Fig 4)

9. Rinse both the brush head and handle under running water wipe them dry. (Fig 5)

10. Put the handle back into the charger.

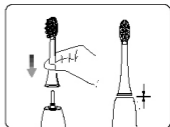
Fig 1

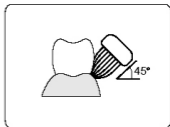
Fig 2

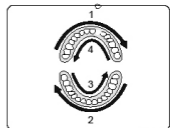
Fig 3

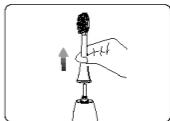
Fig 4

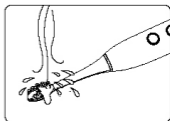
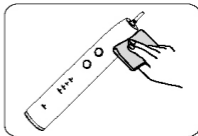
Fig 5


Be sure to remove the plug from the electrical outlet before cleaning the product.

Toothbrush

1. Remove the brush head.
2. Wipe all surface areas with a damp cloth and place the handle back into the charger.
3. Do not submerge the handle into water. Rinse the brush head and bristles after each use.
3. Do not clean in the dishwasher.

**Charger**

1. Using a soft cloth, wipe off any dirt on the charger and the storage.
2. Clean the plug periodically to ensure it remains in good working order. The accumulation of dust or dirt on the plug can cause an electric shock, overheat or fire.

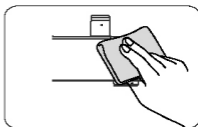**TROUBLESHOOTING**

Problem	Possible cause	Possible cause
The toothbrush doesn't work.	You have just purchased the product or it has not been used for more than 3 months.	Charge the main unit for at least 24 hours.
Weak or diminished vibration.	The toothbrush is pressed hard against your teeth.	Place the brush head gently against your teeth.
The toothbrush operates only for a short period of time even immediately after being charged.	The battery has reached the end of its life cycle.	Contact service center for support.
	Charging time was insufficient.	Charge the main unit for at least 24 hours.

CORRECT DISPOSAL OF THIS PRODUCT

This marking indicates that this product should not be disposed with other household wastes. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

The power source of this product is lithium battery which can be recycled. At the end of the product's life, always remove the battery before you discard the product, and recycle or dispose of the battery according to local waste management requirements. If you have recycling questions, contact your local waste management office.

Removing and recycling the battery

Do not disassemble the product at any time other than when disposing of it. Doing so earlier will destroy the product. Observe basic safety precautions. The appliance must be disconnected from the supply mains when removing the battery.

1. Do not disassemble the device at any time other than when disposing of it.
2. Doing so earlier will destroy the product.
3. The appliance must be disconnected from the power supply when disassembling the device.

Sonic Toothbrush

1. Before removing the battery, be sure to completely drain the unit of power.
2. Use a screwdriver to remove the screw and battery cover. (Fig.6)
3. Hold the handle upside down, press the toothbrush shaft, and remove the two hooks to release the internal components of the handle. (Fig.7)
4. Insert the screwdriver next to the battery connections, and twist to break the connections. Remove the circuit board and pry the battery from the plastic carrier. (Fig.8)
5. Be careful not to short-circuit the negative and positive poles of the battery. Do not connect the appliance to the mains again after you have removed the battery.
6. Recycle or dispose of the battery properly. Ask your local waste management office if they accept the remaining toothbrush components for recycling.

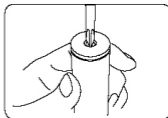
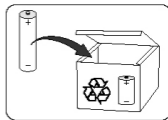
Fig6

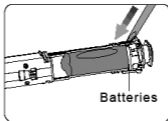
Fig8

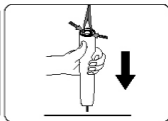
Fig7

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Warranty conditions:

The plant guarantees normal operation of the product within 12 months from the date of the sale, provided that the consumer complies with the rules of operation and care provided by this Manual.

Service life of the product is 5 years. In order to avoid any misunderstanding we kindly ask you to carefully study the Owner's Manual, the terms of warranty, and check the correctness of filling the Warranty Card. The Warranty Card is valid only if the following information is correctly and clearly stated: model, serial number of the product, date of sale, clear seals of the company-seller, signature of the buyer. The model and serial number of the product must correspond to those specified in the Warranty Card.

If these terms are violated, and if the information specified in the Warranty Card is changed, erased or rewritten, the Warranty Card is considered invalid. Setup and Installation (assembly, connection, etc.) of the appliance are described in the attached documentation; they can be carried out both by the user, and by specialists having the necessary qualification sent by USC or selling companies (on a paid basis). In this case, the person (the Institution) installing the appliance shall be responsible for the correctness and quality of installation (setup). Please pay attention to the importance of the correct installation of the appliance, both for its reliable operation, and for obtaining warranty and free service. Require the installation specialist to enter all the necessary information about the installation in the Warranty Card.

In case of failure of the appliance during the warranty period by the fault of the manufacturer, the owner is entitled to free warranty repair upon presentation of the correctly filled Warranty Card together with the appliance in the factory configuration packed in the original packaging to the warranty workshop or the place of its purchase. Satisfaction of the customer's claim due to the fault of the manufacturer shall be made in accordance with the law «On Protection of Consumers' Rights». In the case of warranty repair the warranty period shall be extended for the period of repair and shipment.

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WARRANTY CARD

ARDESTO

Dear Buyer! Congratulations on your purchase of the ARDESTO brand appliance, which was designed and manufactured in accordance with the highest quality standards, and we thank you for choosing this particular instrument.

We ask you to keep the coupon during the warranty period. When purchasing a product, require a full warranty card.

1. Warranty service is carried out only if there is a correctly and clearly filled original warranty card, which indicates: product model, date of sale, serial number, warranty service period, and the seller's seal. *
2. The service life of household appliances is 5 years.
3. The product is intended for consumer use. When using the product in commercial activities, the seller / manufacturer does not bear warranty obligations, after-sales service is performed on a paid basis.
4. Warranty repair is carried out within the period specified in the warranty card for the product in an authorized service center on the conditions and terms determined by applicable law.
5. The product is withdrawn from the warranty in case of violation by the consumer of the operating rules set forth in the instruction manual.
6. The product is removed from warranty service in the following cases:
 - misuse and non-consumer use;
 - mechanical damage;
 - damage caused by the ingress of foreign objects, substances, liquids, insects;
 - damage caused by natural disasters (rain, wind, lightning, etc.), fire, domestic factors (excessive humidity, dust, aggressive environment, etc.)
 - damage caused by non-compliance of power and cable network parameters with state standards and other similar factors;
 - when operating equipment in the power supply network with a missing single ground loop;
 - in case of violation of seals installed on the product;
 - lack of serial number of the device, or inability to identify it.
7. The warranty does not cover consumables and accessories, as well as filters, shelves, drawers, stands and other containers for storing products.
8. The warranty period is 12 months from the date of sale.

* Tear-off maintenance tickets are provided by an authorized service center.

The completeness of the product is checked. I have read the terms of the warranty service, no complaints.

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