REMINGTON | REVEAL

DIAMOND MICRODERMABRASION

GUARANTEE ⁺¹
2
yr
MD3000

Thank you for buying your new Remington® product. Please read these instructions carefully and keep them safe. Remove all packaging before use.

IMPORTANT SAFETY INSTRUCTIONS

WARNING – TO REDUCE THE RISK OF BURNS, ELECTROCUTION, FIRE OR INJURY TO PERSONS:

- 1 This appliance can be used from 18 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children.
- 2 An appliance should never be left unattended when plugged into a power outlet, except when charging.
- 3 Do not use this appliance if it is not working correctly, if it has been dropped or damaged, or dropped in water.
- 4 Do not twist or kink the cable, and don't wrap it around the appliance.
- 5 Do not use attachments other than those that we supply.
- 6 Make sure the power plug and cord do not get wet.
- 7 Do not plug or unplug the appliance with wet hands.
- 8 Keep the charger away from water and heated surfaces.
- 9 Never drop or insert any objects into any opening on the device.
- 10 Do not operate where aerosol spray product is being used or where oxygen is being administered.
- 11 Do not operate under a blanket or pillow. Excessive heating can occur and cause fire, electrocution or injury to persons.
- 12 Do not carry this appliance by the cord or use the cord as a handle.
-  **WARNING:** Do not use this appliance near bathtubs, showers, basins or other vessels containing water.
- 13 This appliance is not intended for commercial or salon use.
- 14 Use only as directed.
- 15 This device is designed for dry use only. Do not use wet or in the shower.
- 16 Always unplug the appliance from the electrical outlet before cleaning. To disconnect, turn the device off, then remove plug from outlet.
- 17 After use, the diamond tips can be removed from the device and rinsed under water. Wipe the handle with a dry cloth only for cleaning.
- 18 This device should be supplied with an approved safety isolation adaptors PA-0510E for Europe and PA-0510U for the United Kingdom with the output capacity of 5.0V DC, 1A (adapter outlet).

WARNINGS AND SAFETY PRECAUTIONS

Before you start using the REVEAL Microdermabrasion device, be sure to read all the warnings and safety information to check the product is suitable for you.

CONSULT A DOCTOR BEFORE USING THIS DEVICE IF:

- You are currently undergoing clinical treatment for a specific skin condition, or have any health concerns.
- If you suffer from circulatory problems, or suffer from sensitive skin.
- If you are susceptible to allergic reactions of the skin.

DO NOT USE THIS DEVICE ON:

- Inflamed, irritated, or sunburned skin.
- Dry, sore, sensitive or chapped skin.
- Any areas of skin with dermatological problems.
- Skin which is broken or infected.
- Skin with spots, pimples, warts, moles or varicose or thread veins.
- Tattoo covered skin or with permanent makeup.

The device is designed to be easy to use and pain free. If you experience a noticeably unacceptable level of discomfort, or skin irritation, stop use immediately and assess you are using the device in accordance with the instructions. If you are satisfied you are using the device correctly and the discomfort continues or an extreme reaction occurs, do not use and consult your doctor.

- Individuals with pacemakers and pregnant women should consult a doctor before using this appliance. Not recommended for use by Diabetics.
- Do not use on an infant, invalid or on a sleeping or unconscious person. Do not use on sensitive skin or on a person with poor blood circulation.
- If the product feels excessively hot, switch off and contact Remington Customer Service.
- Hands must be dry when operating the controls or when handling the plug and electrical wall sockets.
- The REVEAL Microdermabrasion Device is not a professional medical device and is not intended to treat medical conditions of the skin. Do not use this device as a substitute for medical attention.
- This device should not be used by any individual suffering from any ailments that would limit the user's ability to operate the device as intended.

The REVEAL Microdermabrasion Device is NOT recommended for use on any of the following skin conditions.

- Active Rosacea
- Extremely fragile capillaries
- Vascular Lesions
- Acne
- Herpetic Lesions, cold sores
- Warts
- Open wounds
- Psoriasis
- Lupus
- Erythematous
- Diabetes Mellitus
- Skin Cancer
- Sunburned Skin
- Those taking blood thinners


GETTING ACQUAINTED WITH THE REVEAL MICRODERMABRASION DEVICE

WHAT IS MICRODERMABRASION?

- Microdermabrasion is the process of removing the outer most layer of dead skin which in turn stimulates the growth of a fresh new skin layer which restores skin radiance and improves skin tone.
- The Microdermabrasion process utilizes diamond grit tips to loosen the dead skins cells coupled with vacuum suction to pull the dead skin cells off the skin.
- The vacuuming process has a secondary function of stimulating blood flow at the outer skin level.


ABOUT YOUR DEVICE

- The REVEAL Microdermabrasion Device includes three different diamond tips to enable you to personalise your treatment. Fig. 1.

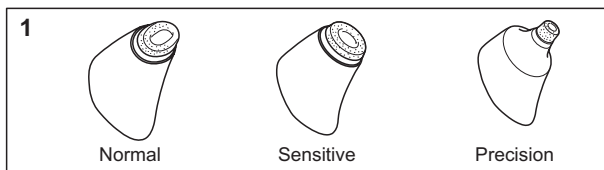- **NORMAL** Diamond Tip – For non-sensitive areas of the face and neck.
- **SENSITIVE** Diamond Tip – For delicate areas or sensitive skin.
- **PRECISION** Diamond Tip – For precise treatment of smaller areas such as the nose and other harder to reach parts.

The REVEAL Microdermabrasion Device also includes four vacuum level settings to allow the user to choose safe, comfortable and effective combinations that work for their skin type.

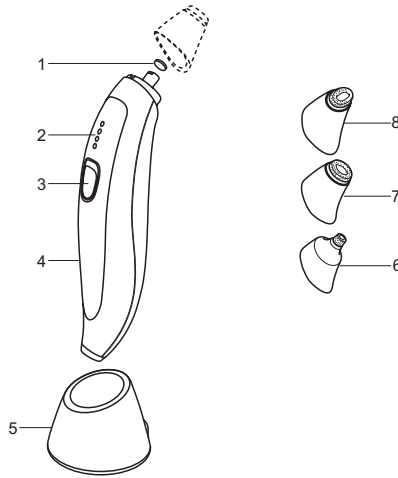KEY FEATURES


- | | |
|------------------------------|------------------------------------|
| 1. Replaceable Filter | 7. Sensitive Diamond Tip |
| 2. LED Power Level Indicator | 8. Normal Diamond Tip |
| 3. Power Level Switch | 9. Replacement Filters (Not Shown) |
| 4. Ergonomic Handle | 10. Storage Pouch (Not Shown) |
| 5. Charging Stand | 11. Power Adaptor (Not Shown) |
| 6. Precision Diamond Tip | 12. Cleaning Brush (Not Shown) |

GETTING STARTED

Be patient when first using your REVEAL Microdermabrasion Device as with any new product, it may take a little while to familiarise yourself with it. Take the time to acquaint yourself with your Microdermabrasion Device, as we are confident you will receive years of enjoyable use and complete satisfaction.

INSTRUCTIONS FOR CHARGING

- 1 Always ensure that your hands, Microdermabrasion Device, Charging Stand and Adaptor are dry before charging.
- 2 Before using for the first time, fully charge the Microdermabrasion Device.
- 3 A full charge will take 4 hours.
- 4 Ensure the device is switched off.
- 5 Insert the Charging Adaptor into the base of the Charging Stand. Fig 2.

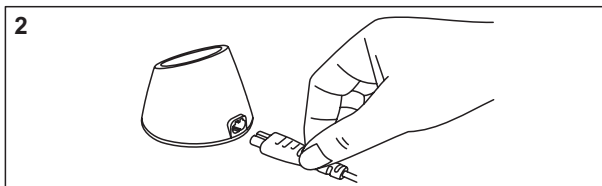- 6 Plug the Charging Adaptor into the wall socket.
- 7 Place the Microdermabrasion Device onto the Charging Stand. The device will need to be oriented so the keyed socket at the end of the device fits into the plug of the Charging Stand.

- 8 The LED lights will begin to light sequentially to indicate the Microdermabrasion Device is charging.
 - 9 When the Microdermabrasion Device is charged, all the LED lights will remain solid for a minute and then will go off to indicate a full charge. A full charge will take 4 hours.
 - 10 One full charge will provide up to 90 minutes of cordless usage time depending on the speed setting selected.
- ✳ NOTE: If the battery charge is low, all four LED lights will flash in unison to indicate the device needs to be recharged.

✳ **QUICK CHARGE FUNCTION:**

- The Microdermabrasion Device is also equipped with a quick charge function.
- To utilise this function, charge the Microdermabrasion Device for 5 minutes.
- This will provide approximately 1 or 2 uses.

◆ INSTRUCTIONS FOR USE

The device is suitable for dry use only. Do not use in the shower or bath.


- ✳ **NOTE: Depending on your skin type, as the Microdermabrasion process utilises diamond tips coupled with vacuum suction, you may experience some reddening or increased sensitivity after your treatment. This is considered perfectly normal and usually diminishes within a few hours.**
- ✳ TIP: As reddening can occur after treatment, it is recommended to carry out the treatments at the end of the day.
- It is therefore strongly recommended that you perform a patch trial on the inside of your forearm to feel how the device works and to see the reaction of the skin.
 - ✳ TIP: It is also advised to carry out this first trial treatment with the Sensitive Diamond Tip using the lowest setting.
 - Once familiar with the Device and its effect on your skin, you can personalise your treatment by adjusting the settings and selecting the appropriate tip.
 - In use you will feel the Diamond Tip against your skin and a gentle suction sensation from the vacuum setting.
 - Adverse reactions are rare, but if skin becomes unduly inflamed or irritated, stop using the device immediately.

◆ HOW TO USE

- 1 Remove all packaging from the device.
- 2 For your convenience, the device is fitted with a filter for your first treatment.
- ✳ TIP: It is recommended that you always check that there is a filter in place before you commence treatment.
- 3 Select your chosen diamond tip to begin your treatment routine.
- ✳ NOTE: The tips are clear wrapped for cleanliness and need to be removed before use. They attach onto the device magnetically and must be fully engaged to provide the best suction performance.
- 4 Wash and towel dry your face before use.
- 5 Pull your hair back from your face and secure with a head band or pony tail.
- 6 Turn the device on by pressing the power switch.
- 7 Select the desired suction level:
 - The Microdermabrasion Device will automatically start on the lowest suction level.
 - The suction levels are adjusted by repeated presses of the power switch.
 - The LED lights will light up to indicate what suction level the device is on. One LED light indicates the lowest setting. Four LED lights indicate the highest setting.
- 8 Press and hold the switch again when you want to turn the device off.

👁️ TREATMENT ROUTINE

- The treatment routine for your face will take about 3 minutes to complete.
- Plan to spend 20 to 30 seconds in each area, following the suggested stroke direction as indicated by the arrows. Fig. 3.

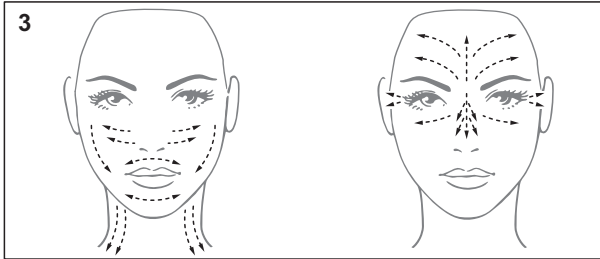- ✳️ NOTE: The exact treatment time will vary depending on your skin type and which zone you are treating.
- **Do not repeatedly move the diamond tip over the area you are treating e.g. after moving the diamond tip over an area once, move to the next untreated area.**
- For best results ensure that the diamond tip of the device is flat and in full contact with the skin.
- Do not press the diamond tip hard against the skin, but allow the diamond tip to do all the work for you.
- During treatment you should feel a gentle suction effect as the diamond tip moves across the skin.
- It is normal to notice the slight coarseness of the diamond tip in use, but it should never cause discomfort.

✳️ TREATMENT TIPS:

- Select the correct Diamond Tip for the area of the face you are treating.
- Keep the diamond suction tip in contact with the skin at all times and ensure to keep the tip moving.
- ✳️ NOTE: Do not allow the tip to stop in the same spot as the prolonged suction could result in pulling blood to the skin surface and result in excessive reddening, soreness or bruising.
- It may help to hold your skin taut with the other hand during treatment.
- Never use the high suction settings on the delicate skin areas around your eyes.
- Never use the device directly underneath your eyes or on your eyelids.
- Do not use the device over moles, tattoos and permanent make-up.
- Always use the diamond suction tip in the recommended direction as shown on the illustration.

✳️ POST TREATMENT:

- Rinse your face with cold water to remove any additional dead skin cells.
- Apply your favourite nourishing moisturizing cream to the entire treated area.
- Avoid using acid based moisturizing creams.
- Apply sunscreen to the treatment area if you plan to go outdoors. Avoid sun bathing soon after treatment.

✳️ FREQUENCY OF TREATMENTS:

- With Microdermabrasion, it is important to allow enough time between treatments for new skin cells to regenerate on the surface of the skin to reveal healthy looking and glowing skin. This will vary from person to person depending on age and skin type.
- It is also important to allow time between treatments so you can assess your skin and determine the appropriate time for your next treatment.
- **For this reason, it is recommended not to exceed one or two treatments a week.**


👁 FITTING AND REPLACING FILTERS

The Microdermabrasion treatment process removes dead skin cells from the outer surface of your skin. The device is equipped with disposable filters that collect the skin particles during treatment.

- For your convenience the device is fitted with a filter for your first treatment and replacement filters are included in the pack.
- For optimum hygiene and product effectiveness, it is important to replace the filter after every use and to place a new and unused filter into the device before each treatment.
- ✳ **NOTE:** Be certain there is always a fresh filter in the device while in use. Use of the device without a filter would draw the skin particles into the pump mechanism and shorten the life of the device and invalidate the warranty.

✳ TO REPLACE THE FILTER:

- Remove the Tip from the Handle.
- Using your fingertip, remove the circular black filter from the end of the suction tube and discard.
- Place a new filter into the suction tube and gently press it into place with your fingertip. Fig. 4.

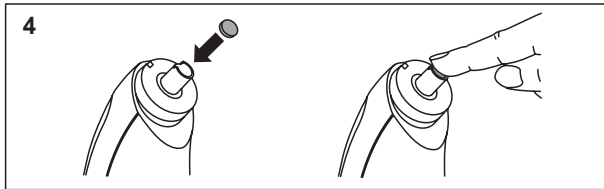
👁 CHANGING THE DIAMOND TIPS

The tips are clear wrapped for cleanliness and this needs to be removed prior to first use. They attach onto the device magnetically and must be fully engaged to provide the best suction performance.

✳ Removing or Changing the Diamond Tips

- Hold the Microdermabrasion Device with one hand
- With the other hand, gently pull the Diamond Tip upwards until it comes off. Fig 5.
- To replace, gently position the Diamond Tip and the magnetic force will pull it into place. Fig 6.

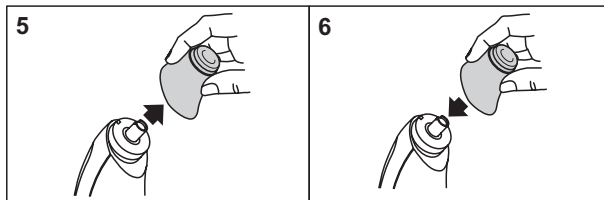👁 CLEANING AND MAINTENANCE

- ✳ **NOTE:** The handle of the Microdermabrasion Device is not waterproof. Do not submerge the handle or hold it under running water to clean.
- To ensure long lasting performance and cleanliness of your Microdermabrasion device, clean the diamond tips and handle after each use.
- Ensure the Microdermabrasion Device is switched off before cleaning.
- Remove the diamond tip and brush away excess skin cells with the cleaning brush provided. When removed, the tips can be rinsed under running water to clean. Allow to air dry.
- ✳ **NOTE:** The Diamond Tips may accumulate skin cells over time, reducing their effectiveness. If this happens, then the Diamond Tips can be periodically cleaned using water enhanced with an alcohol based barber's cleaning solution.

- Use caution with the diamond tips to prevent damage. Avoid dropping the tips. Do not use harsh chemicals, scouring pads or abrasive cleaning agents on the tips.
- Do not clean the device or any of the components in the dishwasher.
- ✳ **CAUTION:** Ensure all parts are dry before connecting the charger into a wall socket.

✳ **STORAGE**

- Store your device in a clean, dry place and out of the reach of children.
- Do not store away while the device is wet from cleaning. Allow to air dry before storing.
- Avoid storing the abrasive tips where other items can come in contact and damage them.

PRODUCT DISPOSAL

The battery must be removed from this Microdermabrasion Device before it can be disposed of. The appliance must be disconnected from the mains supply when removing the battery.

- ✳ **NOTE:** The battery in this device is not replaceable. When the battery has exceeded its life it must be discarded.
- To remove the rechargeable battery:
- Run the Microdermabrasion Device until the battery is completely discharged.
- Remove the diamond tip from the device.
- Remove the 4 round screw caps at the base of the suction tube.
- Remove the 4 screws from the retaining cap.
- Pry the front and rear housing covers apart.
- Remove the square silver battery and power circuit board from the device.
- Dispose of the battery safely.

ENVIRONMENTAL PROTECTION

To avoid environmental and health problems due to hazardous substances in electrical and electronic goods, appliances marked with this symbol must not be disposed of with unsorted municipal waste, but recovered, reused or recycled.

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© 2017 Spectrum Brands, Inc.
Remington® is a Registered Trade Mark of Spectrum Brands, Inc.,
VARTA Consumer Batteries GmbH & Co. KGaA, Alfred-Krupp-Str. 9, 73479 Ellwangen, Germany