ARDESTO™

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HEALTHY TEETH AND GUMS


You can effectively and easily take care of your gums, teeth, and beauty of your smile, using ARDESTO irrigator on a daily basis! With the combination of water pressure and pulsation, the irrigator removes harmful bacteria and food debris deep between teeth and below gumline braces and other orthodontic constructions where brushing and even dental floss won't reach.

Traditional Clean

Traditional brushing can't easily reach hidden food debris and harmful bacteria


Flossing has limited access below the gumline and may damage gums

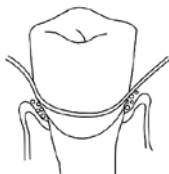IRRIGATOR CLEAN

Removes trapped food debris and harmful bacteria with high frequency pulse water floss.

If you use irrigator on a daily basis, it gives you fresh breath, makes gums and teeth healthier, and helps to clean orthodontic constructions!

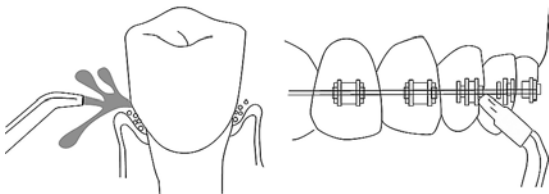Consult with your dentist on the irrigator modes usage.

CONTENTS

IMPORTANT SAFEGUARDS	39
EQUIPMENT	46
TECHNICAL SPECIFICATIONS	46
PRODUCT DESCRIPTION	47
PRODUCT ACCESSORIES	48
OPERATING STEPS	51
USAGE TIPS	58
CLEANING AND MAINTENANCE	59
BATTERY CARE AND DISPOSAL	61
TROUBLESHOOTING	62
DISPOSAL	65

IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, the basic safety precautions below should always be followed.

READ ALL INSTRUCTIONS BEFORE USE DANGERS!

- To reduce the risk of electrocution:
- Do not handle the power cable with wet hands.
- Do not immerse in water or other liquid.
- Do not use while bathing.
- Do not place or store the device where it can fall or be pulled into a tub, shower stall or sink.
- Do not reach for a product that has

fallen into water. Unplug immediately if it was connected to the power supply.

- Do not place the product in a place with high temperature or humidity and avoid direct sunlight. Do not put the product into fire or heat the product.
- Please do not damage, destroy or modify the charging equipment, do not forcefully stretch, bend or twist the charging cable, and do not press heavy objects on the charging cable or let the charging cable get stuck between objects.
- Check the charger cord for damage before the first use and during the life of the product.

To reduce the risk of burns, electric shock, fire, serious injury or damage issues:

- Do not plug this device into a voltage system that is different from the voltage system specified on the device or charger. Do not use this product if it has a broken power cord.
- In case of not working properly, or it has been dropped into any liquid, contact the authorized service center.
- Do not use the device without installing the tip.
- Do not direct water flow under the tongue, into the ear, nose or other delicate organs. This product is capable of producing pressures that may cause serious damage to these organs. Refer to operating

instructions for correct usage.

- Use this product only as indicated in these instructions or suggestions by your dental professional.
- It's normal phenomenon if gum bleed a little when using water flosser for the first time. If bleed much, please stop using and ask for your dentist advice.
- Fill reservoir with water or other solutions recommended by dental professional.
- Do not drop or insert any foreign object into the product.
- Keep the power cord away from heated surfaces.
- Do not operate where oxygen or aerosol sprays are being used.
- Do not use mouth washes/rinses

that contain iodine, salt water or concentrated essential oil insoluble in water in products. Using these liquids will reduce device performance and shorten its service life.

- Oral irrigators shall be used only with cold or warm water, or specific solutions as defined in the manufacturer's instructions
- Do not fill hot water above 40 °C into the water tank to avoid scalding.
- Do not disassemble the product or change the parts of the product to avoid damage or danger.
- Do not use while wearing any oral jewelry. Please put off prior to use.
- Do not use if you have an open wound on your tongue or in your mouth.

- Patients with periodontal disease or those who have had oral surgery in the past 2 months, please consult a doctor before using the device.
- If your physician advise you to receive antibiotic premedication, you should consult your dentist before using this instrument or any other oral hygiene aid.
- For household use only. Do not use outdoors.
- Closely instruct and supervise children and individuals with special needs in the proper use of this product.
- Children should be supervised to ensure that they do not play with the device.
- This device can be used by children

over 8 years of age and persons with physical or sensory or mental disabilities or with insufficient experience and knowledge, if they are under constant supervision or instructed in the safe use of the device and they understand the possible dangers. Cleaning and maintenance by the consumer should not be carried out by children without supervision.


EQUIPMENT

Portative irrigator (with water tank)	1 pc
Water tank	1 pc
Standard tip	2 pcs
Plaque brush tip	1 pc
Orthodontic tip	1 pc
Periodontal tip	1 pc
Tongue cleaner tip	1 pc
USB-A to USB-C cable	1 pc
Travel bag	1 pc
Instruction manual (warranty card included)	1 pc

TECHNICAL SPECIFICATIONS

Rated voltage:	3.7 V DC
Rated power:	8 W
Waterproof rating:	IPX7
Charging time:	≤4 hours
Input parameters:	5V DC, 1A

PRODUCT DESCRIPTION

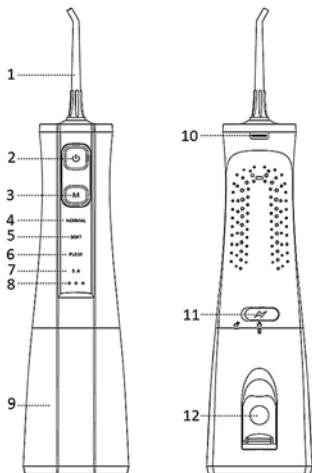- | | |
|--------------------------|----------------------------|
| 1. Tip | 7. SA mode indicator |
| 2. On/Off button | 8. Battery level indicator |
| 3. Mode button | 9. Water tank |
| 4. Normal mode indicator | 10. Tip removal button |
| 5. Soft mode indicator | 11. Charging port |
| 6. Pulse mode indicator | 12. Water tank lid |

PRODUCT ACCESSORIES

Standard tip

Replace every 6 months

- Normal Cleaning

Orthodontic tip

Replace every 6 months

- Braces
- Normal Cleaning

Dental plaque tip

- Dental Implants
- Bridges
- Crowns
- Normal Cleaning

Periodontal tip

- Braces
- Normal Cleaning

Tongue cleaner tip

- Tongue Coating
- Normal Cleaning


Standard tip

Perform a full 360° rinse of hard-to-reach parts of the mouth that can't be cleaned with a toothbrush. More effective than dental floss and healthier for the gums.

Orthodontic tip

The orthodontic tip features small and dense soft bristles that help brush and rinse hard-to-reach parts around the braces.

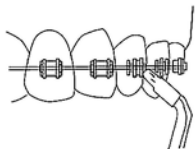Dental plaque tip

The dental plaque tip features three sets of ultra-fine soft bristles that are used to deeply clean hidden areas between the teeth, under the gums, and in cavities. The tip can effectively rinse out debris on the surface of the teeth, bacteria left from food, and even dental plaque.

The device is suitable for users who require special care or deep cleaning, such as those with dental implants, dentures, crowns, or bridges.

Periodontal tip


The periodontal tip features an ultra-fine soft spout and is used to deliver liquid medicine to periodontal pockets.

The device is suitable for users who suffer from gingivitis. Tilt the tip and your teeth to 45°, then gently lower the soft spout to the alveolar. Select the lowest gear, then turn the device on and start cleaning along the gum line. Consult your dentist for detailed instructions on liquid-to-water ratio, dosage, and dilution if necessary.

Tongue cleaner tip

The tongue scraper is used to effectively clean bacteria that coat the tongue surface and solve bad breath. Place the scraper on the back of the tongue coating, select the lowest gear then turn on the device. Scrape the tongue coating while gently working forward.

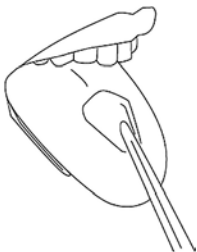
OPERATING STEPS

Very important! Before first use, fully charge the product by connecting the USB cable to the charging port and plugging it into a power source. The product cannot be used while charging.

Step 1. Fill water into water tank

Method 1: open the water tank lid and fill the water tank with water. If the water tank lid comes off, align the gasket of the lid with the tank opening and press it in to close.

Method 2: rotate the water tank clockwise and pull it down to remove it, then fill the water tank with water. After filling the water, align the water tank with the snap-fit point and rotate counterclockwise to reinstall it.

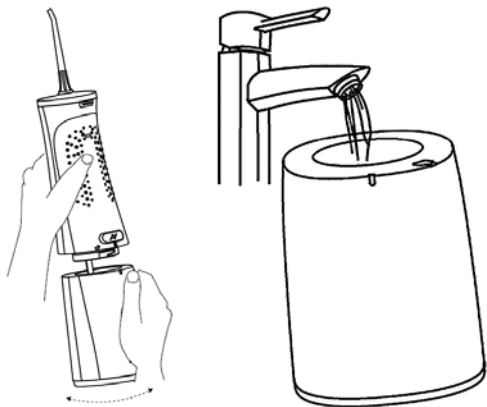
Note: ensure the product is turned off when the water tank is filled with water.

Step 2. Install the tip

Insert the tip into the hole at the top of the main unit. When you hear a 'click' sound, the tip is

installed successfully. Rotate the tip to adjust the direction of the water outflow. To replace the tip, press and hold the tip removal button, and then

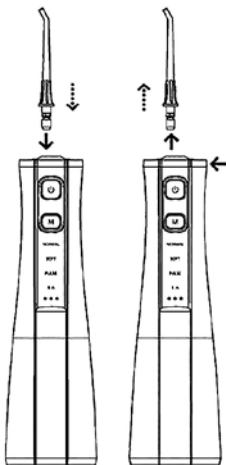take out the tip.

Step 3. Select a mode

This product is available in four modes. When it is turned off, you can select the desired mode by

briefly pressing the mode button. The functions of each mode are as follows:

Press and hold the switch button to activate the working mode. Briefly press the mode button to cycle through the four modes: Normal -> Soft -> Pulse -> SA in sequence. In Normal, Soft and Pulse modes, briefly press the switch button to turn off the device.

SA mode. In SA mode, press and hold the switch button to activate the SA function. Release the switch button to stop the SA function. If the device is not in use, the product will automatically shut down with the indicator light dimming after 5 seconds of standby.

The number of indicator lights illuminated during charging represents the current battery level. The lights will flash in sequence. Once the product is fully charged, all lights will be on. Press and hold the mode button to display the current battery level, indicated by the number of lights on. During use, the indicator light flashes to remind the user to charge. When the indicator light flashes three times, the product will automatically shut down.

Memory function. When the product is turned on, it remembers the mode used before the last shutdown.

Keep the device upright during use to ensure optimal working condition. The device will automatically shut off after running for 120 seconds in any mode.

Step 4. Preparation before use

Before use, lean slightly over the wash basin, insert the tip into your mouth, aim it at your teeth, gums, or slits between teeth, rotate the tip to adjust the water outflow direction, ensuring the water flow is perpendicular to the rinsing area, and keep your mouth slightly open, allowing the rinsing water to flow out of the oral cavity.


Note: when using the oral irrigator for the first time, start with the soft mode. Once your oral cavity adapts, you can use other modes. During use, you may experience itching or slight gum bleeding, which is normal and typically stops within 1 to 2 weeks.

Step 5. Use

Press and hold the switch button for 1 second to start the product. If you need to pause during use, gently press the switch button. To continue using, press and hold the switch button again. This device features 120-second automatic shutdown. The water in the tank will be used up before automatic shutdown. To continue using, gently press the switch button to turn off, add water, and then restart the device. Please turn off the device after use.

Step 6. After use

After use, open the water tank, pour out the remaining liquid in the tank, turn on the product, and discharge the residual liquid inside the product until no more water sprays out, then turn off the

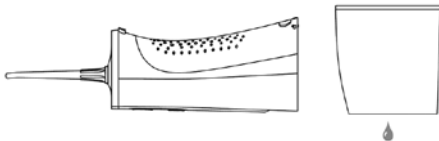device.


Step 7. Charging

The indicator light is white and remains on during

use. When the indicator light flashes 3 times and the product shuts down, it indicates a low battery state. When the product is in a low battery state, pressing the switch button will not activate it. The indicator light will flash 4 times, signaling the need for immediate charging. The indicator light flashes while charging and stays on when the product is fully charged. The indicator light remains on for 20 minutes before off.

USAGE TIPS

For optimal cleaning results, clean towards the front teeth from the cheek tooth (back teeth). When cleaning, move the tip along the gum line at a 90-degree angle and pause briefly between teeth. When using, slightly close your lips to prevent water from splashing, leaving a small gap to allow water in the mouth to flow freely into the wash basin.

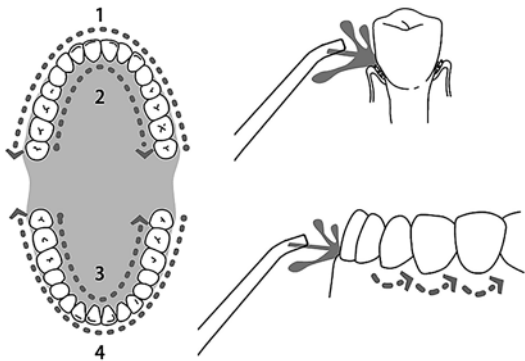CLEANING AND MAINTENANCE

Cleaning the body

1. Rotate the water tank counterclockwise and pull it down to remove the water tank.
2. Rinse the inside of the water tank with room temperature water; do not use any detergent for cleaning.
3. After cleaning, align the water tank with the snap-fit point and rotate clockwise to reinstall it.

Cleaning the water tank

1. Hard water deposits may accumulate on the product, depending on the mineral content in the water. If not cleared in time, they will reduce device performance.

2. Aim the tip at the wash basin, turn on the product, and let it run until the water is used up. Refill the water tank with clean warm water and repeat the above operation for cleaning. Repeat this cleaning operation every 1 to 3 months to ensure device performance.

BATTERY CARE AND DISPOSAL

Battery Maintenance

- The battery must be fully charged before first use.
- For long-term storage, please fully charge the battery before storage.
- To extend the battery lifespan, please charge it immediately when the power is low.

Battery Recycling

This product uses a rechargeable battery that can be recycled. Before discarding the product at the end of its life, be sure to remove the battery and recycle or dispose of it according to local waste management regulations.

When the product is turned off, remove the tip and the bottom cover of the water tank, use a screwdriver or other tools to pry open the bottom cover of the main body and take out the built-in battery for recycling.

Note: do not short-circuit the positive and negative terminals of the removed battery. Apply tape to the lead terminals for insulation.

TROUBLESHOOTING

Problem	Cause	Solution
Leakage between the tip and the main unit	The seal ring inside the tip holder is damaged, or the tip is not properly installed	Install the tip in its place.
Contact the authorized service center	Heavy wear on the soft nozzle	Replace of the periodontal pocket nozzle with a new one
Water pump inlet leakage	Seal ring damaged or impurities at the inlet edge	Contact the authorized service center
Water tank leakage	Seal ring damage Replace	Replace the seal ring. Contact the authorized service center
Water leakage at the connection between the water tank and the main unit due to tilting	The seal ring of the main unit base is deformed, damaged, or missing	Contact the authorized service center to replace the base seal ring
Insufficient pressure	Use in a low battery state	Use after fully charged
The product cannot start	Use in a low battery state or other malfunction	Use after fully charged. Contact the authorized service center

Problem	Cause	Solution
The product indicator light is normal but does not start	Motor jammed or other technical issue	Contact the authorized service center
No water flowing from the product	The product is used for the first time or not used for a long term	Add an appropriate amount of water to the tip installation port, and install the tip before use
	Foreign objects in the tip hole cause blockage	Remove foreign objects or replace the tip
	The filter screen at the connector of the water pumping silicone tube is clogged with foreign objects	Remove the foreign objects
	The water pumping silicone tube is broken	Replace the silicone tube. Contact the authorized service center to replace the silicone tube
Product auto power-on	Water exposure to the control board (water ingress into the mechanism)	Contact the authorized service center

The product cannot be charged/charged properly	The adapter has no output voltage or the charging cable has poor contact	Replace the adapter or charging cable
	The product does not support charging cables with Type-C connectors on both ends (excluding the 'figure-8' connector style)	Use the charging cable provided with the product or an optional charging cable with a USB-A connector
	Charging base dirty/oxidized causing poor contact	Regular cleaning; keep the charging port clean and dry
Automatic shutdown	The product features automatic shutdown	Restart the device
	Battery is depleted	Charge the device
Noise	Noise level ≤ 63 decibels within 100 cycles during product operation with load is considered normal	-

DISPOSAL

The packaging is made entirely of recyclable materials, which you may dispose of at local recycling facilities.

Observe the marking of the packaging materials for waste separation, which are marked with abbreviations (a) and numbers (b) with following meaning:

1–7: plastics / 20–22: paper and fiberboard / 80–98: composite materials.

Dispose of the packaging materials environmentally. The packaging is made entirely of recyclable materials which should be disposed of at local recycling facilities.

To help protect the environment, please dispose of the product properly when it has reached the end of its useful life and not in the household waste.

Information on collection points and their opening hours can be obtained from your local authority.

Contact your local refuse disposal authority for more details of how to dispose of your worn-out product.

The manufacturer reserves the right to change the design and the specifications of the unit without a preliminary notification.

In case of any printing errors or misunderstandings in the content of this manual, the producer reserves the right of final interpretation. Product upgrades or data updates will be included in the new version of the manual without further notice. In case of any change to the product or its appearance color, the actual item inside the package shall prevail.

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This image shows a single sheet of white paper with horizontal ruling lines. The lines are evenly spaced and run across the width of the page. There are no margins, text, or other markings on the paper.

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Warranty conditions:

The plant guarantees normal operation of the product within 12 months from the date of the sale, provided that the consumer complies with the rules of operation and care provided by this Manual.

Service life of the product is 3 years.

In order to avoid any misunderstanding we kindly ask you to carefully study the Owner's Manual, the terms of warranty, and check the correctness of filling the Warranty Card. The Warranty Card is valid only if the following information is correctly and clearly stated: model, serial number of the product, date of sale, clear seals of the company-seller, signature of the buyer. The model and serial number of the product must correspond to those specified in the Warranty Card.

If these terms are violated, and if the information specified in the Warranty Card is changed, erased or rewritten, the Warranty Card is considered invalid.

Setup and installation (assembly, connection, etc.) of the appliance are described in the attached documentation; they can be carried out both by the user, and by specialists having the necessary qualification sent by Authorized Service Center or selling companies (on a paid basis). In this case, the person (the institution) installing the appliance shall be responsible for the correctness and quality of installation (setup). Please pay attention to the importance of the correct installation of the appliance, both for its reliable operation, and for obtaining warranty and free service. Require the installation specialist to enter all the necessary information about the installation in the Warranty Card. This warranty does not cover installation and uninstallation of the appliance. All costs associated with the installation and uninstallation of the appliance (including warranty repairs cases) are borne by the customer.

In case of failure of the appliance during the warranty period by the fault of the manufacturer, the owner is entitled to free warranty repair upon presentation of the correctly filled Warranty Card together with the appliance in the factory configuration packed in the original packaging to the warranty workshop or the place of its purchase. Satisfaction of the customer's claim due to the fault of the manufacturer shall be made in accordance with the law «On Protection of Consumers' Rights». In the case of warranty repair the warranty period shall be extended for the period of repair and shipment.

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WARRANTY CARD

ARDESTO™

Dear Buyer! Congratulations on your purchase of the ARDESTO brand appliance, which was designed and manufactured in accordance with the highest quality standards, and we thank you for choosing this product.

Keep this card during the warranty period. When purchasing a product ask for the warranty card to be completely filled out.

1. Warranty service is carried out only if there is a correctly and clearly filled original warranty card, which indicates: product model, date of sale, serial number, warranty service period, and the seller's seal. *
2. The service life of household appliances is 3 years.
3. The product is intended for consumer use. When using the product in commercial activities, the seller / manufacturer does not bear warranty obligations, after-sales service is performed on a paid basis.
4. Warranty repair is carried out within the period specified in the warranty card for the product in an authorized service center on the conditions and terms determined by applicable law.
5. The product is withdrawn from the warranty in case of violation by the consumer of the operating rules set forth in the instruction manual.
6. The product is removed from warranty service in the following cases:
 - misuse and non-consumer use;
 - mechanical damage;
 - damage caused by the ingress of foreign objects, substances, liquids, insects;
 - damage caused by natural disasters (rain, wind, lightning, etc.), fire, domestic factors (excessive humidity, dust, aggressive environment, etc.)
 - damage caused by non-compliance of power and cable network parameters with state standards and other similar factors;
 - when operating equipment in the power supply network with a missing single ground loop;
 - in case of violation of seals installed on the product;
 - lack of serial number of the device, or inability to identify it.
7. The warranty does not cover consumables and accessories as well as blades, filters, shelves, drawers, stands and other containers for storing products, etc.
8. The warranty period is 12 months from the date of sale.

* Tear-off maintenance tickets are provided by an authorized service center.

The completeness of the product is checked. I have read the terms of the warranty service, no complaints.

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